The comprehensive range of bone punches from Aesculap!



Aesculap Neurosurgery



KERRISON Bone Punches from Aesculap



Non-detachable KERRISON Bone Punch

Proven

...considered a classic for over 100 years. The graceful slim design and excellent cutting properties are its hallmark.

Aesculap – a tradition of high quality!

Aesculap has manufactured KERRISON bone punches for over 100 years.

Throughout the world, "Made by Aesculap" is synonymous with high quality and durability.

Detachable KERRISON Bone Punch

Detachable

...can be disassembled without tools in only 3 steps for efficient cleaning.

From raw material to KERRISON bone punch.


The high quality of the Aesculap KERRISON bone punches is a result of the longstanding experience and quality awareness of our employees.




Skilfully and expertly hand-crafted

In order to meet the high demands of every operation, every single part is manufactured by hand and subsequently checked.

Non-detachable KERRISON punches

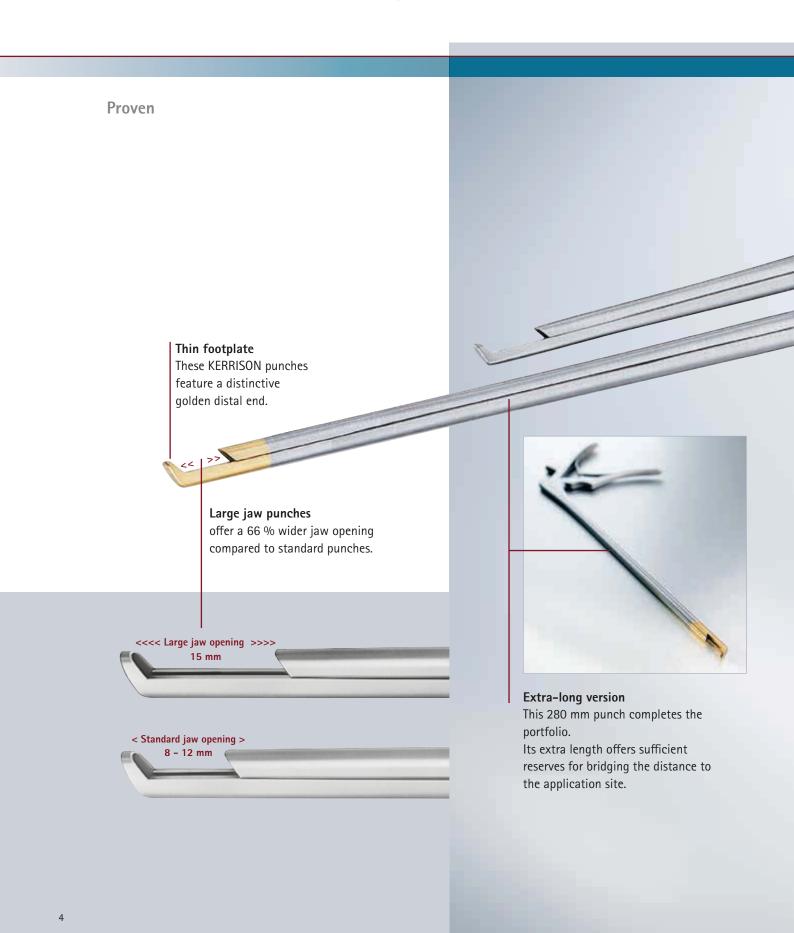

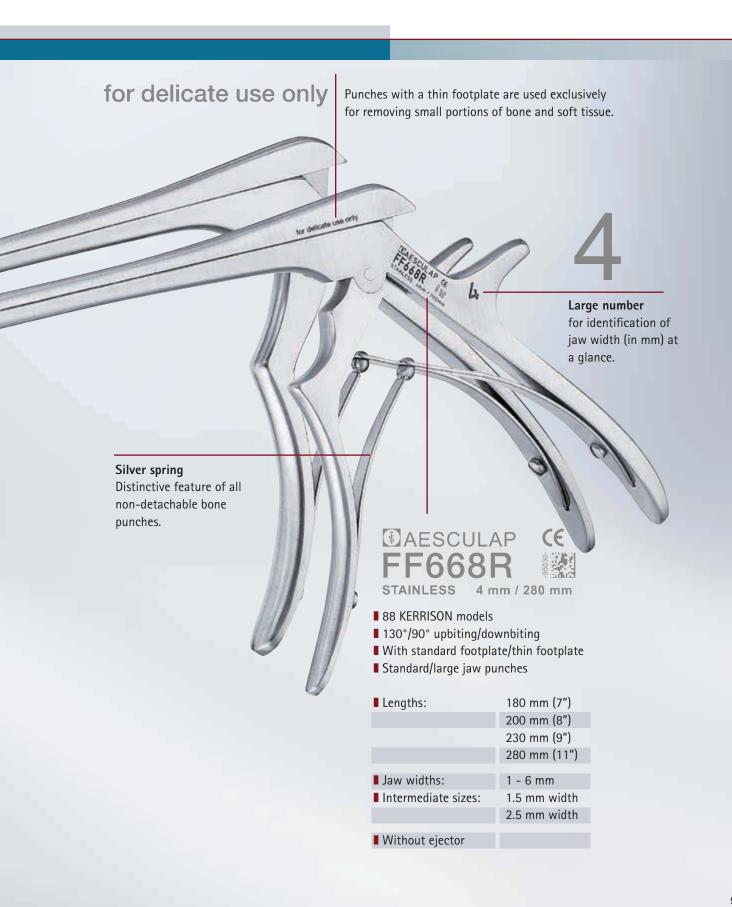


Non-detachable KERRISON punches

Jaw position 130°, upwards							
Footplate							
		Star	ndard	Sta	ndard	TI	hin
				Jaw o	pening		
		Standard	(mm)	Large	(mm)	Standard	(mm)
Length	Width	_					
180 mm	1.0 mm	FF715R	(8)	_	-	FF771R	(8)
	1.5 mm	FF638R	(9)	FF647R	(15)	FF645R	(9)
	2.0 mm	FF717R	(9)	FF677R	(15)	FF772R	(9)
	2.5 mm	FF639R	(10)	FF648R	(15)	FF646R	(10)
	3.0 mm	FF720R	(10)	FF678R	(15)	FF773R	(10)
	4.0 mm	FF724R	(12)	FF679R	(15)	FF769R	(12)
	5.0 mm	FF722R	(12)	FF680R	(15)	FF770R	(12)
	6.0 mm	FF726R	(12)	FF681R	(15)	-	-
200 mm	1.0 mm	FF691R	(8) New	-	-	FF661R	(8)
	1.5 mm	FF649R	(9)	-	-	FF658R	(9)
	2.0 mm	FF737R	(9)	FF660R	(15)	FF662R	(9)
	2.5 mm	FF657R	(10)	-	-	FF659R	(10)
	3.0 mm	FF738R	(10)	FF739R	(15)	FF663R	(10)
	4.0 mm	FF727R	(12)	FF745R	(15)	FF664R	(12)
	5.0 mm	FF728R	(12)	FF751R	(15)	FF665R	(12)
	6.0 mm	FF729R	(12)	FF753R	(15)	-	_
230 mm	1.0 mm	FF674R	(8) New	-	-	FF734R	(8) New
	2.0 mm	FF682R	(9)	-	-	FF746R	(9)
	3.0 mm	FF683R	(10)	_	-	FF747R	(10)
	4.0 mm	FF684R	(12)	_	-	FF748R	(12)
	5.0 mm	FF685R	(12)	_	-	FF749R	(12)
280 mm	1.0 mm	FF696R	(8) New	-	-	FF673R	(8) New
	2.0 mm	FF686R	(9) New	-	-	FF666R	(9)
	3.0 mm	FF687R	(10) New	-	-	FF667R	(10)
	4.0 mm	FF688R	(12) New	-	-	FF668R	(12)
	5.0 mm	FF689R	(12) New		-	FF669R	(12)

Jaw position 1	30°.	downwards)
----------------	------	-----------	---

		Footplate					
		Sta	Standard		Standard		hin
				Jaw	opening		
		Standard	(mm)	Large	(mm)	Standard	(mm)
Length	Width						
180 mm	1.0 mm	FF711R	(8)	-	-	FF781R	(8)
	2.0 mm	FF712R	(9)	-	-	FF782R	(9)
	3.0 mm	FF713R	(10)	-	-	FF783R	(10)
	4.0 mm	FF714R	(12)	-	-	FF784R	(12)
	5.0 mm	FF725R	(12)	-	-	FF785R	(12)

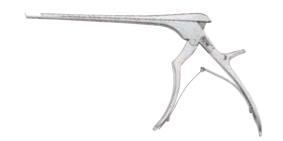

Jaw position 130°, downwards

		Footplate	
		Standard	
		Jaw opening	
		Standard	(mm)
Length	Width		
230 mm	2.0 mm	FF692R	(9)
	3.0 mm	FF693R	(10)
	4.0 mm	FF694R	(12)
	5.0 mm	FF695R	(12)

Jaw position 90°, upwards

		Standard	
		Jaw opening	
		Standard	(mm)
Length	Width		
180 mm	1.0 mm	FF716R	(8)
	2.0 mm	FF718R	(9)
	3.0 mm	FF750R	(10)
	4.0 mm	FF752R	(12)
	5.0 mm	FF754R	(12)
	6.0 mm	FF755R	(12)

Footplate

Jaw position 90°, downwards

		Footplate	
		Standard	
		Jaw opening	
		Standard	
Length	Width		
180 mm	1.0 mm	FF741R	(8)
	2.0 mm	FF719R	(9)
	3.0 mm	FF721R	(10)
	4.0 mm	FF742R	(12)
	5.0 mm	FF723R	(12)
	6.0 mm	FF743R	(12)

All Non-detachable KERRISON punches without Ejector

Width (mm)


1.0 0 1.5 0 2.0 0 2.5 0 3.0 4.0 5.0 6.0 (

Detachable KERRISON punches


Detachable

- Easy disassembly and assembly without tools
- Efficient cleaning and sterilization

Disassembly in 3 steps:

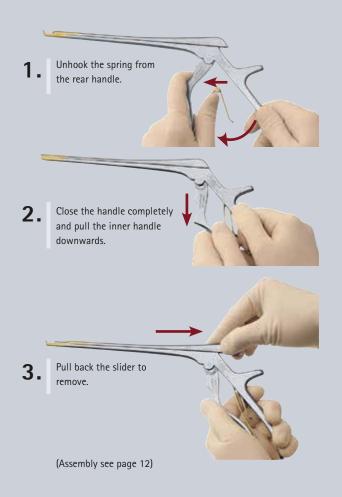



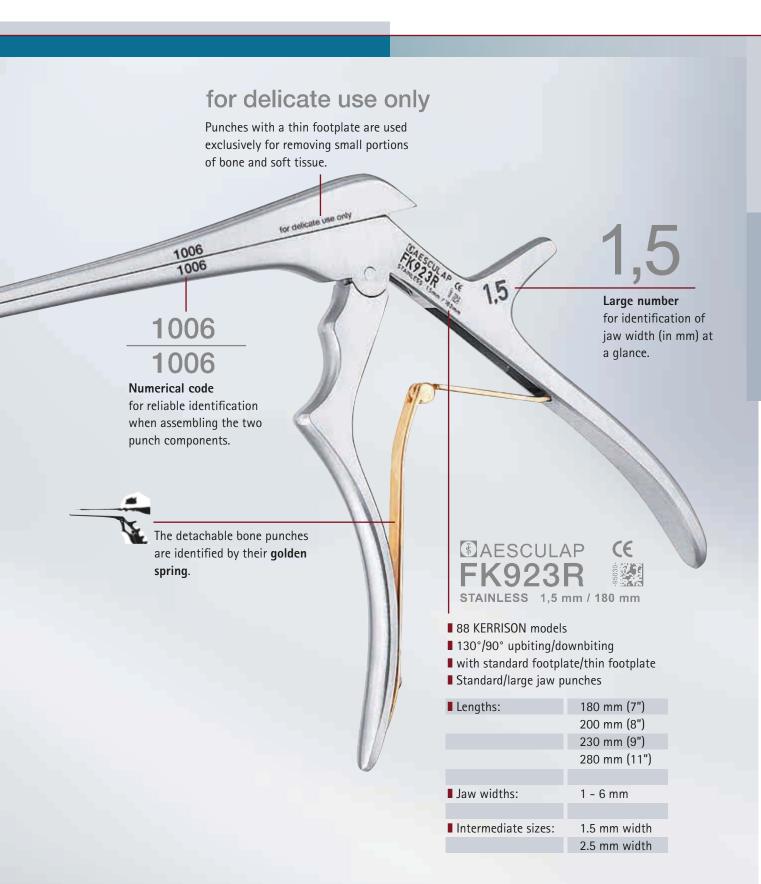

Detachable KERRISON punches

Jaw position 130°, upwards							
Footplate							
		Star	ndard	Sta	ndard	T	hin
				Jaw o	pening		
		Standard	(mm)	Large	(mm)	Standard	(mm)
Length	Width			<u></u>		_	
180 mm	1.0 mm	FK900R	(8)	-	-	FK906R	(8)
	1.5 mm	FK911R	(9)	FK941R	(15)	FK923R	(9)
	2.0 mm	FK901R	(9)	FK971R	(15)	FK907R	(9)
	2.5 mm	FK912R	(10)	FK942R	(15)	FK924R	(10)
	3.0 mm	FK902R	(10)	FK972R	(15)	FK908R	(10)
	4.0 mm	FK903R	(12)	FK973R	(15)	FK909R	(12)
	5.0 mm	FK904R	(12)	FK974R	(15)	FK910R	(12)
	6.0 mm	FK905R	(12)	FK975R	(15)	-	-
200 mm	1.0 mm	FK960R	(8) New	-	-	FK961R	(8)
	1.5 mm	FK966R	(9)	-	-	FK968R	(9)
	2.0 mm	FK913R	(9)	FK970R	(15)	FK962R	(9)
	2.5 mm	FK967R	(10)	-	-	FK969R	(10)
	3.0 mm	FK914R	(10)	FK982R	(15)	FK963R	(10)
	4.0 mm	FK915R	(12)	FK983R	(15)	FK964R	(12)
	5.0 mm	FK916R	(12)	FK984R	(15)	FK965R	(12)
	6.0 mm	FK917R	(12)	FK985R	(15)	-	-
230 mm	1.0 mm	FK918R	(8) New	-	-	FK929R	(8) New
	2.0 mm	FK919R	(9)	-	-	FK925R	(9)
	3.0 mm	FK920R	(10)	-	_	FK926R	(10)
	4.0 mm	FK921R	(12)	-	_	FK927R	(12)
	5.0 mm	FK922R	(12)	-	-	FK928R	(12)
280 mm	1.0 mm	FK981R	(8) New	-	-	FK980R	(8) New
	2.0 mm	FK986R	(9) New	-	-	FK976R	(9)
	3.0 mm	FK987R	(10) New	-	-	FK977R	(10)
	4.0 mm	FK988R	(12) New	-	-	FK978R	(12)
	5.0 mm	FK989R	(12) New	-	-	FK979R	(12)

law	position	130°	. d	lownwards	

		Footplate					
		Standard		Standard		Thin	
			Jaw opening				
		Standard	(mm)	Large	(mm)	Standard	(mm)
Length	Width						
180 mm	1.0 mm	FK930R	(8)	-	-	FK936R	(8)
	2.0 mm	FK931R	(9)	-	-	FK937R	(9)
	3.0 mm	FK932R	(10)	-	-	FK938R	(10)
	4.0 mm	FK933R	(12)	-	-	FK939R	(12)
	5.0 mm	FK934R	(12)	_	-	FK940R	(12)



Jaw position 130°, downwards

		Footplate	
		Standard	
		Jaw opening	
		Standard	(mm)
Length	Width		
230 mm	2.0 mm	FK943R	(9)
	3.0 mm	FK944R	(10)
	4.0 mm	FK945R	(12)
	5.0 mm	FK946R	(12)

Jaw position 90°, upwards

		Footplate	
		Standard	
		Jaw opening	
		Standard	(mm)
Length	Width		
180 mm	1.0 mm	FK948R	(8)
	2.0 mm	FK949R	(9)
	3.0 mm	FK950R	(10)
	4.0 mm	FK951R	(12)
	5.0 mm	FK952R	(12)
	6.0 mm	FK953R	(12)

Jaw position 90°, downwards

		Toutplate	
		Standard	
		Jaw opening	
		Standard	(mm)
Length	Width	—	
180 mm	1.0 mm	FK954R	(8)
	2.0 mm	FK955R	(9)
	3.0 mm	FK956R	(10)
	4.0 mm	FK957R	(12)
	5.0 mm	FK958R	(12)
	6.0 mm	FK959R	(12)


Ejector	×	X	✓	✓	~	~	~	~
Width (mm)	1.0	1.5	2.0	2.5	3.0	4.0	5.0	6.0


Noir® - coated KERRISON punches

Aesthetic

- Easy disassembly and assembly without tools
- Efficient cleaning and sterilization
- Black Noir® coating no irritating light reflections

Assembly in 3 steps:

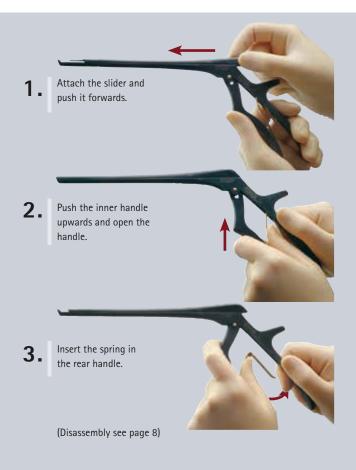

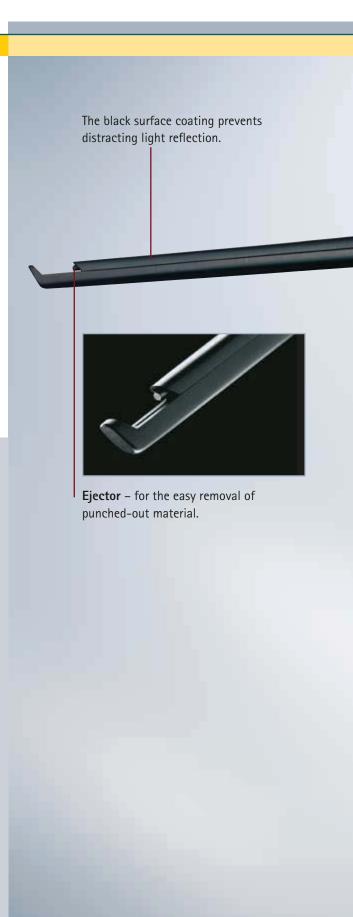




Noir® – coated KERRISON punches

Jaw position 130°, upwards							
	Footplate						
		Star	ndard	Sta	ndard	Thin	
				Jaw o	pening		
		Standard	(mm)	Large	(mm)	Standard	(mm)
Length	Width	_				_	
180 mm	1.0 mm	FK900B	(8)	-	-	FK906B	(8)
	1.5 mm	FK911B	(9)	FK941B	(15)	FK923B	(9)
	2.0 mm	FK901B	(9)	FK971B	(15)	FK907B	(9)
	2.5 mm	FK912B	(10)	FK942B	(15)	FK924B	(10)
	3.0 mm	FK902B	(10)	FK972B	(15)	FK908B	(10)
	4.0 mm	FK903B	(12)	FK973B	(15)	FK909B	(12)
	5.0 mm	FK904B	(12)	FK974B	(15)	FK910B	(12)
	6.0 mm	FK905B	(12)	FK975B	(15)	-	-
200 mm	1.0 mm	FK960B	(8) New	-	_	FK961B	(8)
	1.5 mm	FK966B	(9)	-	-	FK968B	(9)
	2.0 mm	FK913B	(9)	FK970B	(15)	FK962B	(9)
	2.5 mm	FK967B	(10)	-	_	FK969B	(10)
	3.0 mm	FK914B	(10)	FK982B	(15)	FK963B	(10)
	4.0 mm	FK915B	(12)	FK983B	(15)	FK964B	(12)
	5.0 mm	FK916B	(12)	FK984B	(15)	FK965B	(12)
	6.0 mm	FK917B	(12)	FK985B	(15)	-	-
230 mm	1.0 mm	FK918B	(8) New	-	-	FK929B	(8) New
	2.0 mm	FK919B	(9) New	-	-	FK925B	(9)
	3.0 mm	FK920B	(10) New	-	-	FK926B	(10)
	4.0 mm	FK921B	(12) New	-	_	FK927B	(12)
	5.0 mm	FK922B	(12) New	-	-	FK928B	(12)
280 mm	1.0 mm	FK981B	(8) New	-	_	FK980B	(8) New
	2.0 mm	FK986B	(9) New	-	-	FK976B	(9)
	3.0 mm	FK987B	(10) New	-	-	FK977B	(10)
	4.0 mm	FK988B	(12) New	-	-	FK978B	(12)
	5.0 mm	FK989B	(12) New		-	FK979B	(12)

Jaw	position	90°,	upwards
-----	----------	------	---------

		Footplate					
		Standard		Standard		T	hin
		Jaw opening					
		Standard	(mm)	Large		Standard	(mm)
Length	Width						
180 mm	1.0 mm	FK948B	(8)	-	-	-	-
	2.0 mm	FK949B	(9)		-		-
	3.0 mm	FK950B	(10)	-	-	-	-
	4.0 mm	FK951B	(12)	-	-	-	-
	5.0 mm	FK952B	(12)	-	-	_	-
	6.0 mm	FK953B	(12)	_	-		-



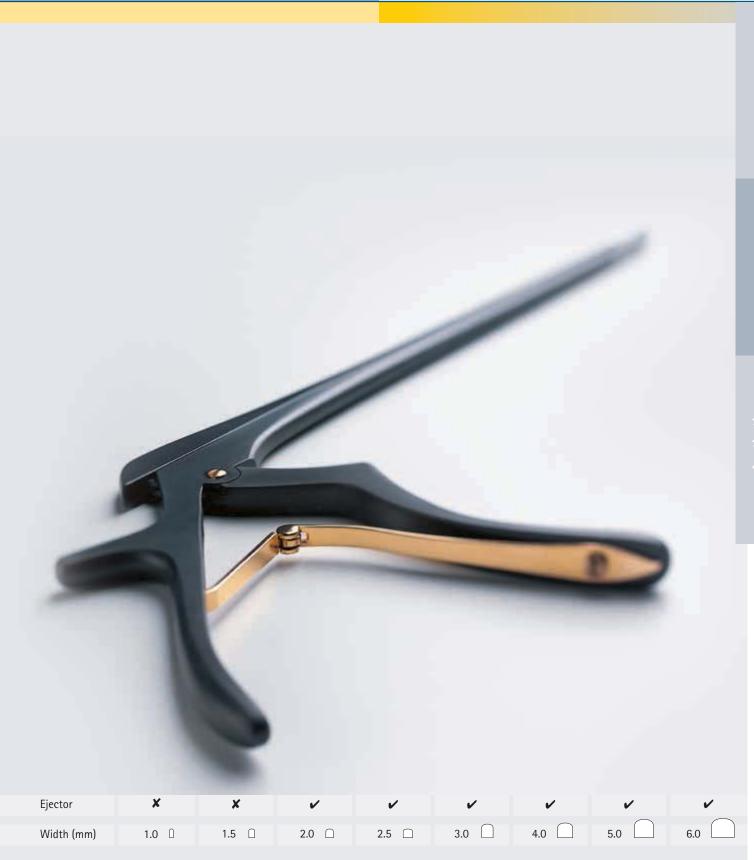

Special KERRISON Bone Punches

Special

In addition to standard punches, Aesculap offers a portfolio of special bone punches for almost every surgical situation.

- 8 special bayonet-shaped and curved punches
 - Classic non-detachable bone punch with an upbiting curved jaw for easy removal of bone material, especially from between the vertebral arches.

Curved KERRISON punch, upbiting

Length	Width	Jaw opening	Art. no.
200 mm	2.0 mm	10 mm	FF466R
	3.0 mm	10 mm	FF467R
	4.0 mm	10 mm	FF468R
	5.0 mm	10 mm	FF469R

• Classic non-detachable bayonet-shaped bone punch for improved vision over the shaft of the instrument in the operating area.

Bayonet-shaped KERRISON punch, 130°, upbiting

Length	Width	Working length	Jaw opening	Art. no.
240 mm	2.0 mm	170 mm	10 mm	FF496R
	3.0 mm	170 mm	10 mm	FF497R
	4.0 mm	170 mm	10 mm	FF498R
	5.0 mm	170 mm	10 mm	FF499R

■ Extra long – 12 KERRISON models




■ Intermediate sizes – 30 KERRISON models



■ Large jaw – 36 KERRISON models





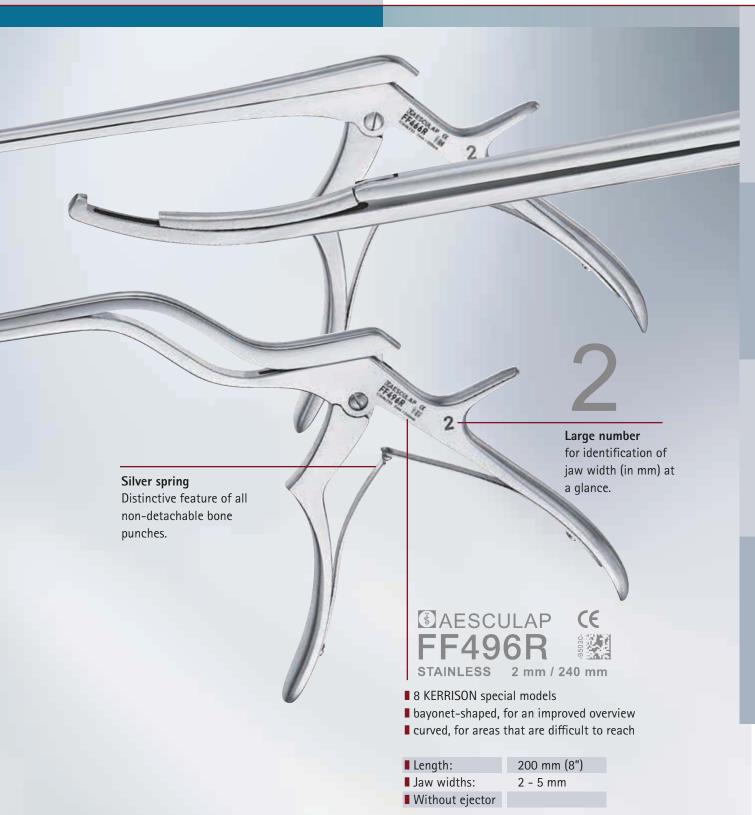


Aesculap® Instrument Organisation System


Good reasons for choosing Aesculap instrument baskets:

- Special Storage frame JF120R for bone punches available
- With large perforation specially suited for cleaning of instruments contaminated with rough debris (e.g. bone chips)
- Smooth surfaces for preventing damage to surgical gloves
- Additional user comfort thanks to low weight construction
- Comprehensive product range offering solutions not only for using instrument baskets in the Sterile Container System but also for single-use packaging materials
- Fully compatible with the Aesculap Instrument Organisation System





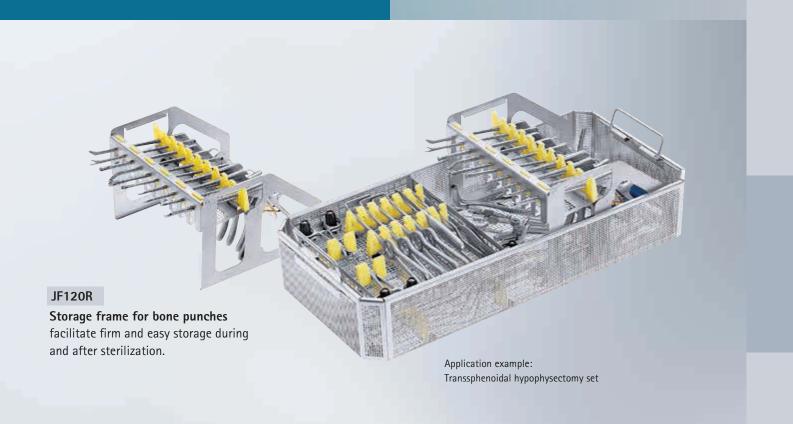

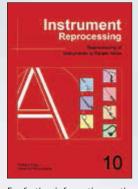

reprocessing of instruments, please refer to the brochure "Proper Maintenance of Instruments".



For further information on the Information on disassembly and reprocessing of the KERRISON bone punch is also available on the KERRISON poster C22911.



For further information on the Aesculap Instrument Organization System, please refer to the brochure C81902.

KAIRison® Pneumatic Bone Punch

The benefits

- Easy handling
- Reduces physical effort required
- Relieves strain on hands and joints
- Less risk of reduced motor skills through fatigue
- Innovative safety mechanism
- Short learning curve
- Excellent balance
- Suitable for right- or left-handed use
- Reliable
- Ergonomic design

The feedback

"My first reaction after using the KAIRison® was to wonder why this unique tool had not been invented 20 years sooner."

Dr. Joseph Maroon, UPMC, Pittsburgh, PA, USA

Technology that makes life easier.





KAIRison® Pneumatic Bone Punch

In just three easy steps



Connect the tube to the compressed air outlet in the OR. Then connect the tube to the pressure reducer, and connect the disposable tube to the handle.



2. Insert the shaft, and turn this toward the handle until you hear it click into place. Please check the shaft is positioned correctly.



3. Switch the KAIRison® on. It is now ready to use. It is that simple.



C25511 KAIRison® short guide poster

Easy shaft replacement



Pull the hock (1) and turn the shaft (2).



7 To detach the shaft just lift it up.





Insert the new shaft (1) and turn this toward the handle until you hear it click into place (2) – "click".



C25411 KAIRison® assembly/ disassembly poster



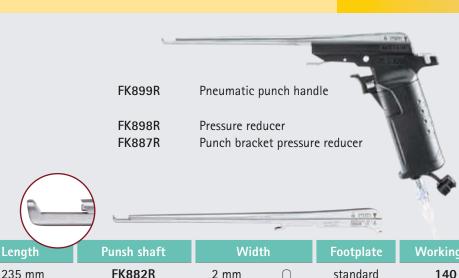

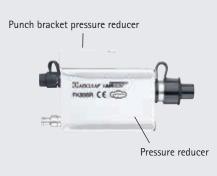

Length	Punsh shaft	Width	Footplate	Working length	Jaw opening	Ejector	
235 mm	FK882R	2 mm 🗆	standard	140 mm	12 mm	✓	
	FK883R	3 mm	standard	140 mm	12 mm	~	
	FK884R	4 mm	standard	140 mm	12 mm	~	
	FK885R	5 mm	standard	140 mm	12 mm	✓	
	FK886R	6 mm	standard	140 mm	12 mm	V	
285 mm	FK892R	2 mm 🗆	standard	190 mm	12 mm	~	
	FK893R	3 mm	standard	190 mm	12 mm	~	
	FK894R	4 mm	standard	190 mm	12 mm	V	
	FK895R	5 mm	standard	190 mm	12 mm	✓	
	FK896R	6 mm	standard	190 mm	12 mm	· ·	

FK888SU


Disposable tube (3 m)

Sales unit: PAK = Package of 10 tubes



FK881R

Storage rack for KAIRison® handle and pressure reducer

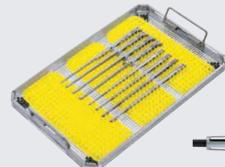


FK891R

Storage rack for KAIRison® shafts



GA466

Pneumatic hose (5 m) with Aesculap-Dräger connection

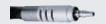

GA468

Pneumatic hose (5 m) with Aesculap-Schrader connection



The main product trademark 'Aesculap' and the product trademarks 'Noir' and 'KAIRison' are registered trademarks of Aesculap AG.

Subject to technical changes. All rights reserved. This brochure may only be used for the exclusive purpose of obtaining information about our products. Reproduction in any form partial or otherwise is not permitted.